

PIN

EN

OVERLAY TF

SIZING CHART

Compatible with distal locking systems, cushion liners (non-locking), and lanyard suspensions.


To determine the correct size, begin by identifying the suspension type, then measure the limb's length and circumference at a specific height.

Smart Size Guide


[Click here](#)


STEP 1

Which suspension ?

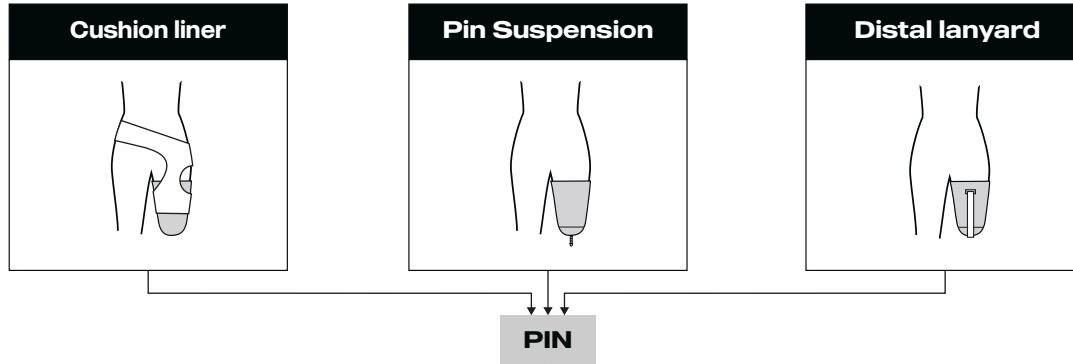

STEP 2

Which length ?

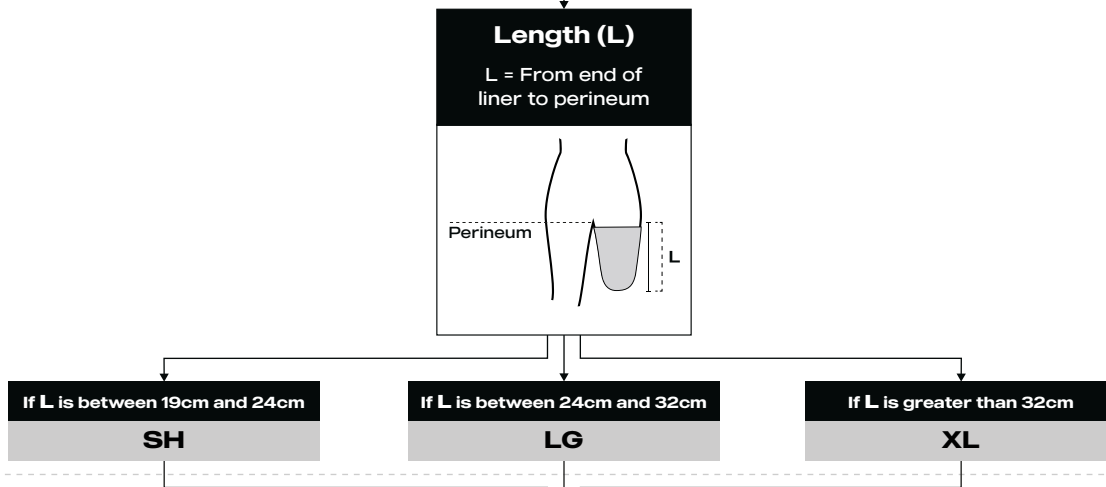


STEP 3

Which size ?

Circumference (C)

Height (H) = 8 cm from distal end

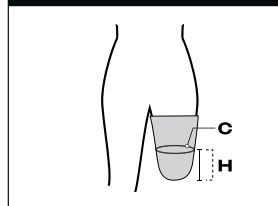


Refer to the chart for the product code based on measured circumference at **8 cm from distal end**.

Suspension	Length (L)	Circumference (C)	Code
Distal Locking	19-24 cm	32-36 cm	OVTF-32-SH-PIN
		36-40 cm	OVTF-38-SH-PIN
		40-44.5 cm	OVTF-40-SH-PIN
		44.5-48 cm	OVTF-44-SH-PIN
		48-52 cm	OVTF-48-SH-PIN
		52-62 cm	OVTF-52-SH-PIN
Distal lanyard	24-32 cm	32-36 cm	OVTF-32-LG-PIN
		36-40 cm	OVTF-38-LG-PIN
		40-44.5 cm	OVTF-40-LG-PIN
		44.5-48 cm	OVTF-44-LG-PIN
		48-52 cm	OVTF-48-LG-PIN
		52-62 cm	OVTF-52-LG-PIN
Cushion liner			

Suspension	Length (L)	Circumference (C)	Code
Distal Locking	+32 cm	32-36 cm	OVTF-32-XL-PIN
		36-40 cm	OVTF-38-XL-PIN
40-44.5 cm		OVTF-40-XL-PIN	
Distal lanyard		44.5-48 cm	OVTF-44-XL-PIN
Cushion liner		48-52 cm	OVTF-48-XL-PIN
		52-62 cm	OVTF-52-XL-PIN

Please see next page for
Vacuum suspensions →

ETHNOCARE™

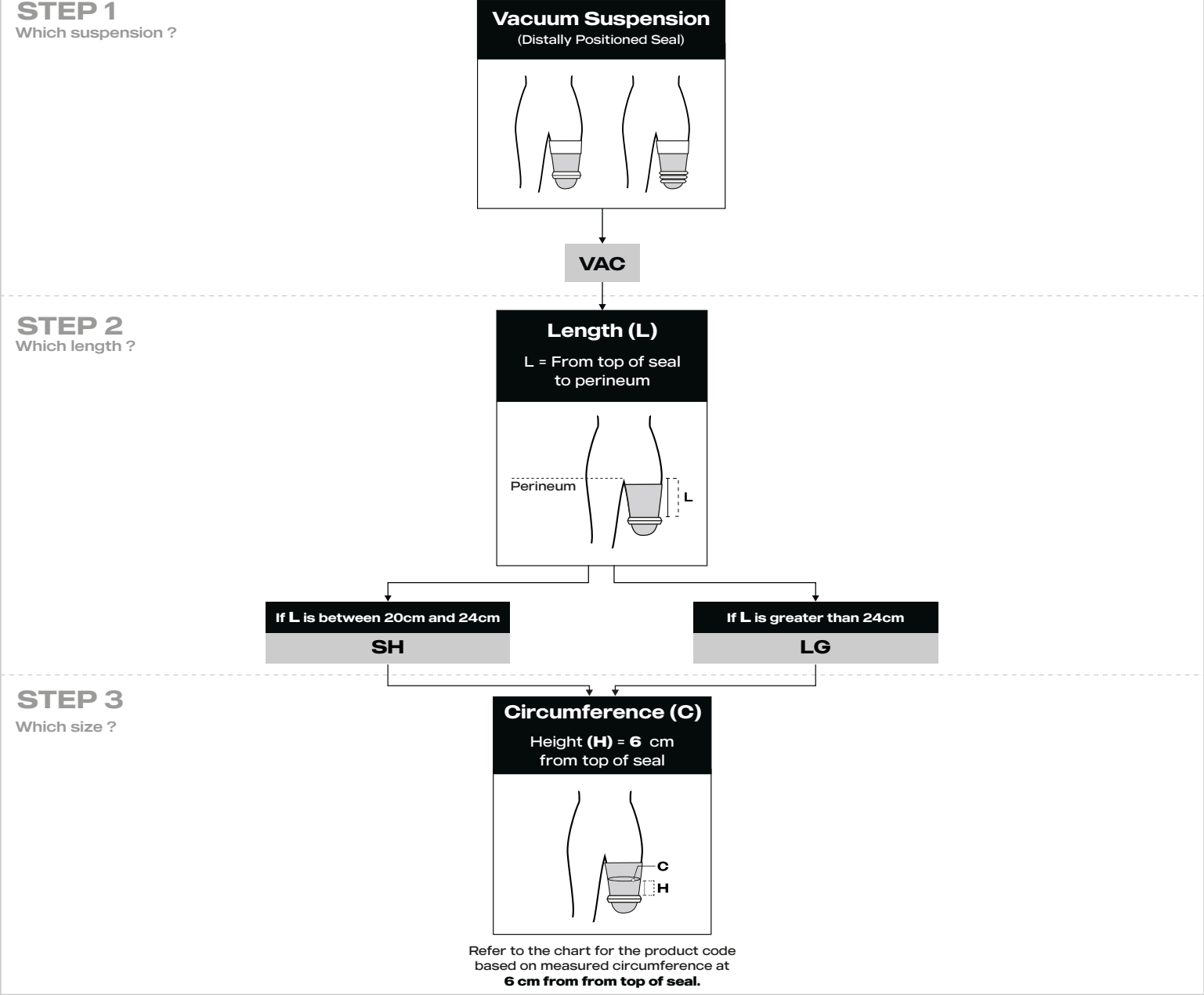

Suspension	Length (L)	Circumference (C)	Code
Vacuum Suspension (Distally Positioned Seal)	20-24cm	32-38 cm	OVTF-32-SH-VAC
		38-42 cm	OVTF-38-SH-VAC
		42-47 cm	OVTF-40-SH-VAC
		47-50 cm	OVTF-44-SH-VAC
		50-54 cm	OVTF-48-SH-VAC
		54-62 cm	OVTF-52-SH-VAC
	+24 cm	32-38 cm	OVTF-32-LG-VAC
		38-42 cm	OVTF-38-LG-VAC
		42-47 cm	OVTF-40-LG-VAC
		47-50cm	OVTF-44-LG-VAC
		50-54 cm	OVTF-48-LG-VAC
		54-62 cm	OVTF-52-LG-VAC

Please see next page for
Other suspensions →